
# Product instructions

## Specifications

Material: Dupont PBT/ABS/TPE	Power Supply: 1xAAA battery
Frequency: 22,000 strokes per minute	
Product size: 25x25x174 mm	Product Standard: GB4706.59-2008
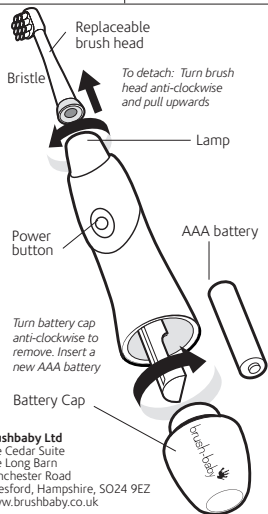
## ENGLISH:

Product contains batteries and/or recyclable electric waste. For environmental protection do not dispose in household waste. For recycling take to electric waste collection points provided by your country.

**Warnings: This product is not a toy. Adult supervision is required at all times. Store out of reach of children. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children. Adult assembly is required. Keep small parts out of reach of children. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been**

**given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.**

Contains 1 x AAA battery. Do not soak the product in water for a long period of time

Do not open, repair or alter the product without authorization. Do not place the product in high temperatures or direct sunlight. The battery is to be inserted with the correct polarity. Non-rechargeable batteries are not to be recharged.

A rechargeable battery is to be removed from the appliance before being charged. An exhausted battery is to be removed from the appliance and safely disposed of.

The battery supply terminals are not to be short-circuited. Do not use attachments which are not recommended by the manufacturer. Replace brush (head) every 3 months

or sooner if brush head becomes worn. For tooth brushing only. Do not chew or bite on brush.

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